

BRAVO™

Chamber Autoclaves



The Perfect Blend of Speed and Capacity

BRAVO is a class B fractionated vacuum sterilizer that incorporates both pre-vac (to remove air) and post-vac (to ensure a dry load every time) cycles, for efficient operation.



Fast

BRAVO 17V and 21V are two of the fastest chamber autoclaves on the market, sterilizing and drying instruments in as little as 24 - 28 minutes.



Gentle

The fresh water feature guarantees fresh water with every cycle. Fresh water reduces the build up of oils and debris.



Versatile

BRAVO offers a number of filling & draining options, including a top and front manual fill and rear automatic fill option that can pull water from a 10L remote bottle placed up to 6 feet away.



Dry

With closed door drying and vacuum technology, the BRAVO is able to dry faster than other chamber autoclaves.



Capacity

BRAVO 17V can sterilize 3 full-size cassettes or up to 20 pouches.

BRAVO 21V can sterilize 3 oral surgery cassettes or up to 25 pouches.



Cycle Times and Maximum Load Weights

Cycles (134°C @ 4 minutes)	Cycle Times (Including Drying**)		Maximum Load Weight	
	17V	21V	17V	21V
Porous (Handpieces) / Wrapped	35 mins.	40 mins.	*	*
Hollow (Handpieces) / Unwrapped	27 mins.	32 mins.	6.0 kg / 13.6 lbs.	7.5 kg / 16.5 lbs.
Solid / Wrapped	24 mins.	28 mins.	3.0 kg / 6.6 lbs.	4.0 kg / 8.8 lbs.
Solid / Unwrapped	19 mins.	23 mins.	6.0 kg / 13.6 lbs.	7.5 kg / 16.5 lbs.

* 17(V) / 21V: Unpackaged porous materials:1kg / 1.25kg, Porous material in single package: 0.75kg / 1.0kg, Solid material/hollow instruments in single pack: 3kg / 4kg, Solid material / hollow instruments in double package:1.5kg / 2.0kg.

** Note: Cycle times represent 220V models

Unit Specifications

	BRAVO 17V	BRAVO 21V
Unit size (L x W x H)†	22.0 x 19.0 x 16.5" / 560 x 480 x 420mm	25.0 x 19.0 x 16.5" / 635 x 480 x 420mm
Volume	4.5 Gal / 17 L	5.5 Gal / 21 L
Chamber Dimensions (D x L)	10 x 13.5" / 254 x 343mm	10 x 17.5" / 254 x 445mm
Weight	128 lbs. / 55 kg	139 lbs. / 63 kg
Voltage	220V, 10A	220V, 10A

†When installed into cabinets, a minimum 4" / 2" (102mm / 51mm) free cool air space is required on the back / sides, and top of the unit.

Recommended Accessories



Ensure™ Monitoring for steam sterilization

An easy to use, in-office monitoring system for steam sterilization.


SCI-BISK-120 - Ensure BI Starter Kit
(25 BI's, Incubator & Log Book)



VistaPure™ AutoFill

Replenish any new or older-model STATIM® and BRAVO autoclave with distilled-quality water automatically, enabling truly hands-free convenience.

S3020 - AutoFill for STATIM
S3030 - AutoFill for BRAVO

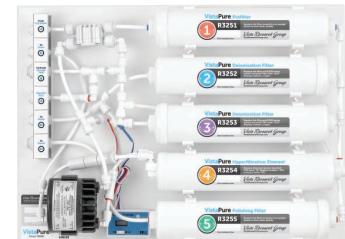

VistaPure Water Filtration

Provides two grades of high-quality water for autoclaves, dental bottles, ultrasonic cleaners and the rinse cycle of instrument washers.

V3000



Ensure Biological Indicators (BI)

Ensure BIs are reliable and self-contained indicators for results in as little as 24 hours. BIs are the only method to ensure effective performance of your sterilizer. The Ensure BIs are available in a starter kit system that consist of 25 BIs, the unique incubator and an easy to fill-in log book.

SCI-BI100 - Biological Indicators (100/box)

Manufactured by:
SciCan Ltd.
1440 Don Mills Rd., Toronto, ON,
Canada, M3B 3P9
T +1 416-445-1600
TF +1 800-667-7733
customerservice@scican.com

VistaPure™
are manufactured by:
Crosstex International, Inc.
1733 W. Parkside Lane
Phoenix, AZ 85027 / USA
T +1 800-531-2221
VistaResearchGroup.com

www.scican.com

COLTENE